

An economical system that allows the accuracy & consistency of advanced technology

State of the art electronics to accurately and efficiently read the Stanbio Uri-Chek[®] 10SG urine strips.



Minimal training time

Data storage

Up to 120 Tests per hour

Lockout functions



STANBIO

LABORATORY

An EKF Diagnostics Company

Quantitative results are better results.
No need to guess at color matching.

Uri-Trak® 120 Urine Analyzer

Includes Analyzer , Strip Holder, Printer
Paper Rolls (2), Fuses, (2.0A)(2), Power
Cord, Quick Guide, User's Guide

Ref: Contents
G1200-001 1 Each

Uri-Chek® 120 10SG Urinalysis Reagent Strips

The Uri-Chek® 120 10SG urinalysis reagent
strip can be read visually or automatically
on the Uri-Trak® 120 to enable point-of-care
clinical decisions.

Ref: Content:
1280-100 100 tests

Data Transfer Cable


Ref: Contents:
G1240-001 1 Each

Uri-Trak® Liquid Urine Controls

Ref: Contents:
1295-210 Urine Control,
Level 1 (10 mL)
Urine Control,
Level 2 (10 mL)

Other accessories:

- Bar code reader G1250-001
- Thermal printer paper G1260-001
- Sticker printer paper G1270-001

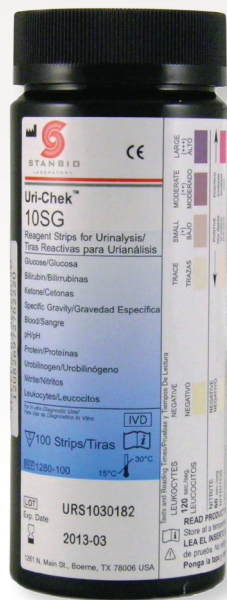

Uri-Trak® 120

- **Minimal training time**
Ease-of-use and user friendly features make testing simple
- **Data storage**
Recalls the last 2,000 test results
- **Lockout Functions**
QC and user lockout programable
- **Built-in thermal printer**
Automatic or manual, up to 3 copies per test, uses thermal or sticker paper
- **Barcode reader available**
Provides quick easy entry of patient data
- **Methodology**
Reflectance photometry
- **Detection**
Photosensitive diode
- **Throughput**
60 tests/hr (single test mode)
120 tests/hr (continuous mode)
- **Strip incubation**
1 min
- **Detection wavelengths**
525 nm and 635 nm
- **Power Source**
100-240V AC, 50/60Hz, 35W
- **Weight**
5.73 lbs / 2.6 kg
- **Dimensions**
(LxWxH) - 10.7" x 10.6" x 5.7"
27cm x 27cm x 14cm

Easy as one,
two, three



Stanbio Laboratory

1261 North Main Street
Boerne, Texas
78006 USA
Tel: 866.782.6246
(U.S. Toll Free)
Tel: 830.249.0772
Fax: 830.249.0851
stanbio@stanbio.com

ISO 13485:2003 Certified



STANBIO
LABORATORY

An EKF Diagnostics Company

RSS0034.01 Eff. Date: 10/09/14 DAR #: 14.0962 ©2014 Stanbio Laboratory